


How to use “FootScan” Foot Scanner

The Scanning of a client’s foot is done in “Fourth parts”.

1. First part is to take image of the foot,
2. Second part is to fill client’s data,
3. Third part is to upload data of client to our website,
4. Fourth part View & print Scanning Report

Part I (Scanning Report)

Start by connecting the USB of the scanner with your Desktop /Laptop computer.
Double Click FootScan icon on desktop. The below Screen will open in which a small window will appear “click connect button”, click **OK**

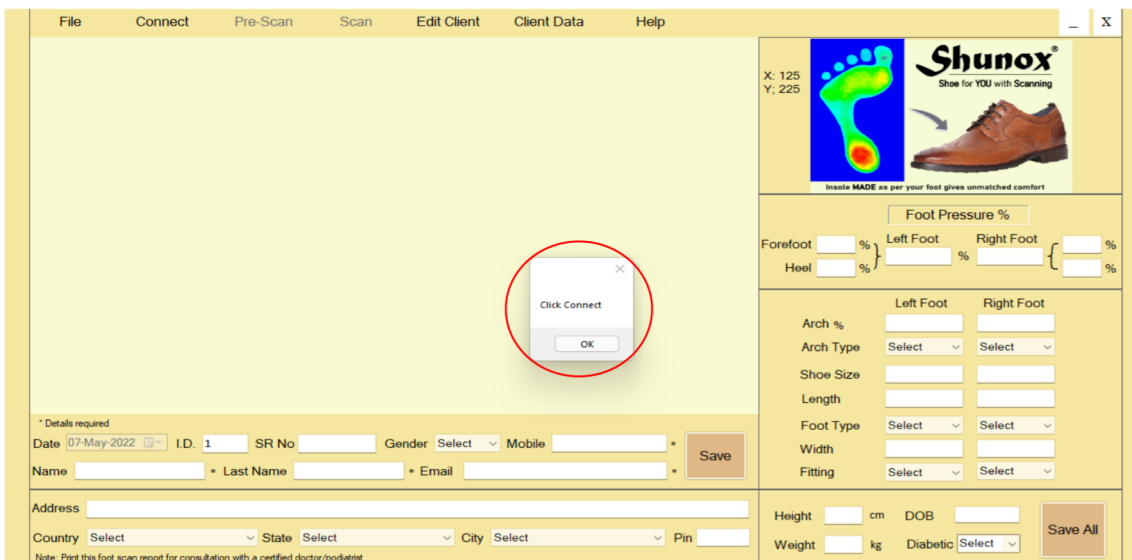


The top left “**Connect**” button will become enabled by turning “Pink”. Click this button.



A small window opens


Please wait...

After a while a small window appears “click Pre- scan button”. Click **OK**

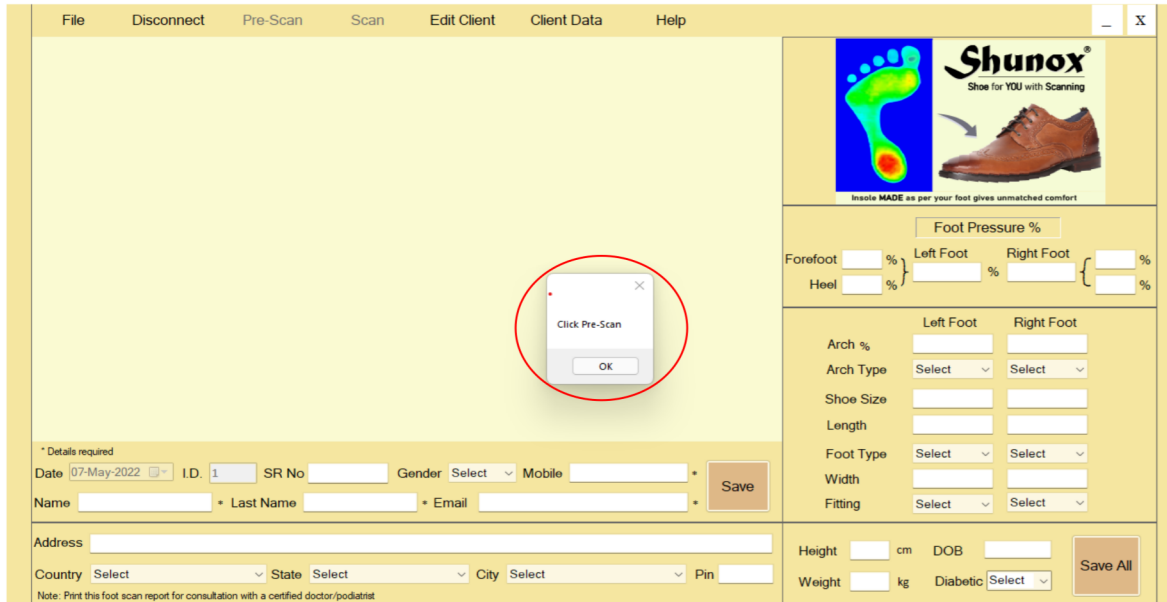

The “Pre-scan” button will become enabled by turning “Pink”. Click this button & a message appears

Please wait...

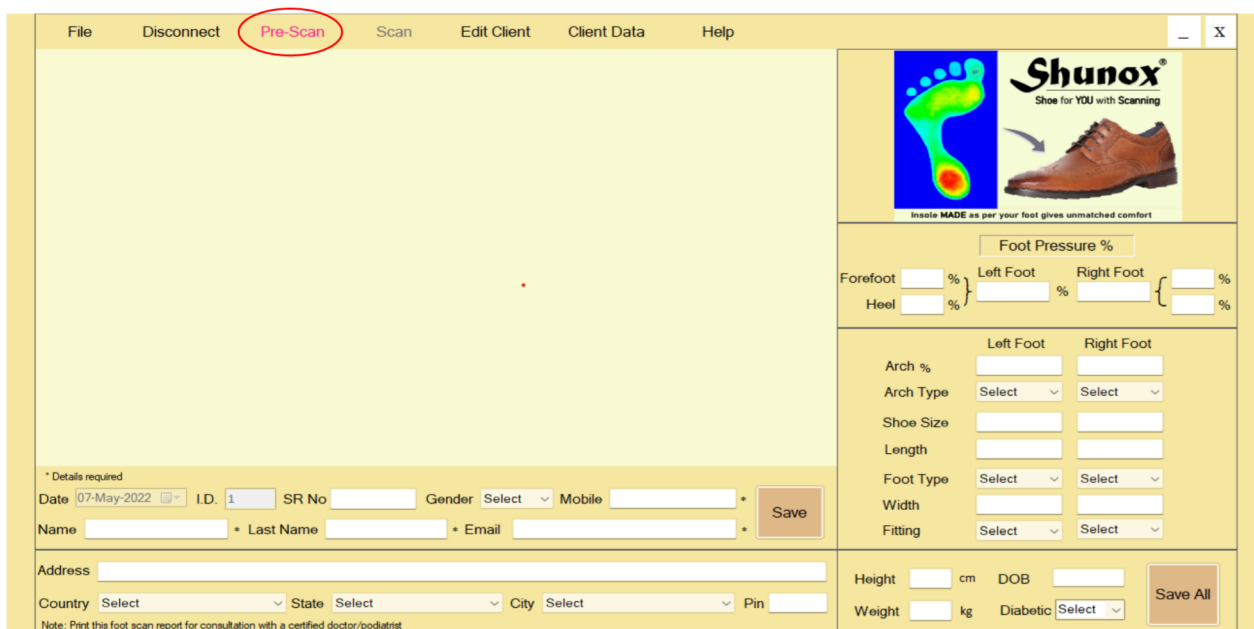

After a moment a small window opens “Please Stand bare foot on the machine”. ASK YOUR CLIENT TO STAND BARE FOOT ON THE SCANNER in an Upright posture & Click **OK**

File Disconnect Pre-Scan Scan Edit Client Client Data Help

Stand bare foot above Scanner.
OK

Inssole **MADE** as per your foot gives unmatched comfort

Foot Pressure %

Forefoot <input type="text"/> %	Left Foot <input type="text"/> %	Right Foot <input type="text"/> %
Heel <input type="text"/> %	<input type="text"/> %	<input type="text"/> %

	Left Foot	Right Foot
Arch %	<input type="text"/>	<input type="text"/>
Arch Type	Select ▾	Select ▾
Shoe Size	<input type="text"/>	<input type="text"/>
Length	<input type="text"/>	<input type="text"/>
Foot Type	Select ▾	Select ▾
Width	<input type="text"/>	<input type="text"/>
Fitting	Select ▾	Select ▾

Height cm DOB

Weight kg Diabetic

Save All

* Details required

Date I.D. SR No Gender Mobile *

Name * Last Name * Email *

Save

Address

Country State City Pin

Note: Print this foot scan report for consultation with a certified doctor/podiatrist

A small window appears “Click scan button”, Click OK

File Disconnect Pre-Scan Scan Edit Client Client Data Help

Inssole **MADE** as per your foot gives unmatched comfort

Foot Pressure %

Forefoot <input type="text"/> %	Left Foot <input type="text"/> %	Right Foot <input type="text"/> %
Heel <input type="text"/> %	<input type="text"/> %	<input type="text"/> %

	Left Foot	Right Foot
Arch %	<input type="text"/>	<input type="text"/>
Arch Type	Select ▾	Select ▾
Shoe Size	<input type="text"/>	<input type="text"/>
Length	<input type="text"/>	<input type="text"/>
Foot Type	Select ▾	Select ▾
Width	<input type="text"/>	<input type="text"/>
Fitting	Select ▾	Select ▾

Height cm DOB

Weight kg Diabetic

Save All

* Details required

Date I.D. SR No Gender Mobile *

Name * Last Name * Email *

Save


Address

Country State City Pin

Note: Print this foot scan report for consultation with a certified doctor/podiatrist

On the top the “**Scan**” button will become enabled by turning “Pink”. Click this button & in a few seconds a “Foot Image of the client” appears on the screen along with foot details & a small window opens “scanning complete”.

File Disconnect Pre-Scan **Scan** Edit Client Client Data Help



Insole **MADE** as per your foot gives unmatched comfort

Foot Pressure %

Forefoot	%	Left Foot	Right Foot	%
Heel	%			

	Left Foot	Right Foot
Arch %		
Arch Type	Select	Select
Shoe Size		
Length		
Foot Type	Select	Select
Width		
Fitting	Select	Select

* Details required

Date 07-May-2022 I.D. 1 SR No Gender Select Mobile * Save


Name * Last Name * Email * Save

Address

Country Select State Select City Select Pin Save All

Note: Print this foot scan report for consultation with a certified doctor/podiatrist

File Disconnect Pre-Scan Scan Edit Client Client Data Help



Insole **MADE** as per your foot gives unmatched comfort

Foot Pressure %

Forefoot	56.62 %	Left Foot	Right Foot	58.36 %
Heel	43.38 %	54.63 %	45.37 %	41.64 %

	Left Foot	Right Foot
Arch %	28.11	17.93
Arch Type	Low	High
Shoe Size		
Length		
Foot Type	Select	Select
Width		
Fitting	Select	Select

* Details required

Date 07-May-2022 I.D. 1 SR No Gender Select Mobile * Save

Name * Last Name * Email * Save

Address

Country Select State Select City Select Pin Save All

Note: Print this foot scan report for consultation with a certified doctor/podiatrist

Confirm

Scanning Completed.
Click [Yes] to proceed further.
Click [No] to "Re-Scan".

Yes No

If you or your client are satisfied with the result, click **YES** & proceed to fill client DATA in the boxes given below.

“BUT” if you or your client feels that the image appearing on the screen is not correct, then ask your client to remain standing on the scanner & click **NO** button.

The top “**Re-scan**” will become enabled by becoming “Pink”. Click this button to again scan foot of your client once again & follow earlier process.

A small window appears stating Data Saved successfully. Click **OK** to complete this foot scanning.

The screenshot shows the Shunox software interface. On the left, a large blue area displays a foot scan result with a color-coded pressure map. A small dialog box in the center reads "Data Saved Successfully" with an "OK" button circled in red. On the right, there's a section for "Foot Pressure %" with data for Left and Right feet. Below this are fields for Arch %, Arch Type, Shoe Size, Length, Foot Type, Width, and Fitting. At the bottom, there's a "Save All" button circled in red. The top menu bar includes File, Disconnect, Pre-Scan, Scan, Edit Client, Client Data, and Help. The bottom section contains client details like Date, I.D., SR No, Gender, Mobile, Name, Last Name, Email, Address, Country, State, City, and Pin.

Foot Pressure %	
Forefoot 56.62 %	Left Foot 54.63 %
Heel 43.38 %	Right Foot 45.37 %
	58.36 %
	41.64 %

	Left Foot	Right Foot
Arch %	28.11	17.93
Arch Type	Low	High
Shoe Size	41	41
Length	258.3	258.3
Foot Type	M+	M+
Width	95	95
Fitting	Medium	Medium

Client Details:

Date: 07-May-2022, I.D.: 1, SR No: Arjun, Gender: Male, Mobile: 9015808124, Name: Sanjay, Last Name: Rana, Email: ranasanjay111@gmail.com, Address: E-189, Sector - 63, Country: India, State: Uttar Pradesh, City: Noida, Pin: 201301.

Part III (Data uploading automatically)

After click Save All Button data uploading are start & please wait sometime. A small window appears stating “Data synching is completed”. Click **OK** to complete this foot scanning

The screenshot shows the Shunox software interface after clicking the "Save All" button. A dialog box in the center reads "Data Synching is completed." with an "OK" button circled in red. The right side of the interface shows the same "Foot Pressure %" and client details sections as before, but the "Save All" button is now disabled. The top menu bar and bottom client details section remain the same.

Foot Pressure %	
Forefoot %	Left Foot %
Heel %	Right Foot %
	%
	%

	Left Foot	Right Foot
Arch %		
Arch Type	Select	Select
Shoe Size		
Length		
Foot Type	Select	Select
Width		
Fitting	Select	Select

Client Details:

Date: 07-May-2022, I.D.: 1, SR No: , Gender: Select, Mobile: , Name: , Last Name: , Email: , Address: , Country: Select, State: Select, City: Select, Pin: .

Part IV (Client Data View or Print)

Click “client data” button on top right side of the screen to open “Upload scanning” page. In this page you will find all foot scanning reports of people. Here you can Print/Delete details

of all clients in local machine & successfully uploaded scanning reports will shift lower side the center line.

[illegible]

This completes the entire scanning process & you can close the application now.